



**Wagner
Grade
1 & 2 DFUs**

143
Wounds

10
Sites


A Composite Ovine Forestomach Matrix and Hyaluronic Acid CAMP for the Treatment of Full Thickness Wounds of the Foot in People with Diabetes: A Randomized Controlled Trial

STUDY PURPOSE:

To evaluate the safety and efficacy of Symphony™ plus standard of care compared with SOC alone in the treatment of Wagner Grade 1 and 2 diabetic foot ulcers.

KEY FINDINGS


57% More DFUs Achieved Complete Closure at 12 Weeks

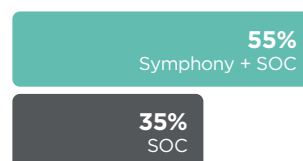

Complete Wound Closure at 12 Weeks

Adjusted PP population: $p=0.039^*$, $n=109$
A 57% relative increase

7.8 Day Reduction in Mean Time to Heal

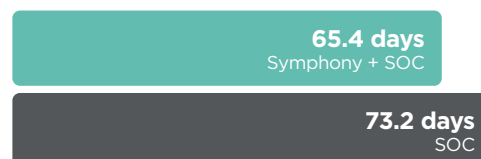

Mean Time to Complete Wound Closure

PP population: $p=0.041^*$, $n=109$
A difference of 7.8 Days



David G. Armstrong, Dennis P. Orgill, Robert D. Galiano, John C. Lantis, Alexander M. Reyzelman, Robert J. Snyder, Adam L. Isaac, Marissa J. Carter, Charles M. Zelen". International Wound Journal. 2026 August; Volume X (Issue Y)

APPLICATION OF SYMPHONY IN WAGNER 2 DFU

Results may vary



Initial Defect

Patient: 63-Year-Old, Male

Wound Location:
Plantar Medial Forefoot

Wound Measurement:
2.2 cm x 1.7 cm x 0.2 cm

Treatment Plan: Weekly sharp debridement and application of Symphony with offloading boot



Week 1

Symphony integrated
Significant reduction in wound area and depth
No complications



Week 4

Symphony integrated
No complications



Week 7 - HEALED

Fully epithelialized
Seven consecutive weekly applications of Symphony

APPLICATION OF SYMPHONY IN WAGNER 1 DFU

Results may vary



Initial Defect

Patient: 58-Year-Old, Male

Wound Location:
Hallux

Wound Measurement:
1.1 cm x 0.8 cm



Week 1

Symphony integrating, anything firmly attached to the wound was left intact and unadhered product was gently debrided
No complications



Week 2

Symphony integrated
Significant reduction in wound area noted
No complications



Week 4 - HEALED

Fully epithelialized
Three consecutive weekly applications of Symphony



Defining the Next Standard in Outpatient Wound Care.
Click or scan the QR code to learn more about Symphony™

aroa.com

For product questions, sampling needs, or detailed clinical questions concerning our products in the US, please call 1-877-627-6224 or email customerservice@aroa.com

RX Only. Prior to use, be sure to read the entire Instructions For Use package insert supplied with the product, or online at <https://eifu.aroa.com/>. AROA™, Aroa Biosurgery™, Symphony™ and AROA ECM™ are trademarks of Aroa Biosurgery Limited. All other trademarks are properties of their respective owners. MKT.2453.00 | ©Sep 2026

